DR.BEI Sonic Electric Toothbrush

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Read the user manual carefully before using and
keep it with care for future reference.

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Contents

| | |
|------------------------------|----|
| Important to note | 1 |
| I. Introduction | 3 |
| II. How to use | 5 |
| III. Specifications | 8 |
| IV. Troubleshooting | 9 |
| V. After-sales service | 10 |
| VI. Environmental protection | 11 |

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Important to note

Statement

First time use may accompany itch or slight bleeding due to the fact that your gum just starts getting used to the high frequency cleaning of electric toothbrush. Kindly suggest begin with Gum-care Mode or Sensitive Brush Head. Itch or bleeding should stop within 1-2 weeks. If not, please discontinue the use and consult your dental professional. This product is not suitable for users that under the age of 6.

★ Tips

1. Don't share brush head with others.
2. Don't surround the charging cable around the base while charging.
3. Don't immerse the travel case into water or other liquid in case of intake.

★ Notes

1. This product is designed for oral cleansing. Do not use it for any other purpose. If you experience any discomfort or pain, please stop using the product and consult your dentist.
2. Check to see if any obvious damage before using, stop using if it damaged or severe fall.
3. Check if the voltage matches with that product specified before charging.
4. If charging base damage, stop using it.
5. Replace brush head every 3 months or sooner if brush head become worn.


★ Warnings

1. Consult your dental professional if you have had orthodontic treatment, dental prosthesis or implantation.
2. This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer prior to use.
3. This product contains a built-in rechargeable battery. Please do not expose the product to flames, and avoid charging, using, or placing it in a high-temperature environment.

4. Don't immerse the charging base in water or other liquid.
5. When charging, make sure the power plug and interface are completely dry before connecting to the power supply.
6. Don't knock or drop the product, stop using if damaged in crash or drop.
7. Don't alter, dismantle or repair the product personally.
8. Don't use it while shower or bath.
9. Don't place metal objects on the base while charging.


I. Introduction

Thank you for choosing DR.BEI to assist on your oral health. Please see the list below to check if you have all the accessories.

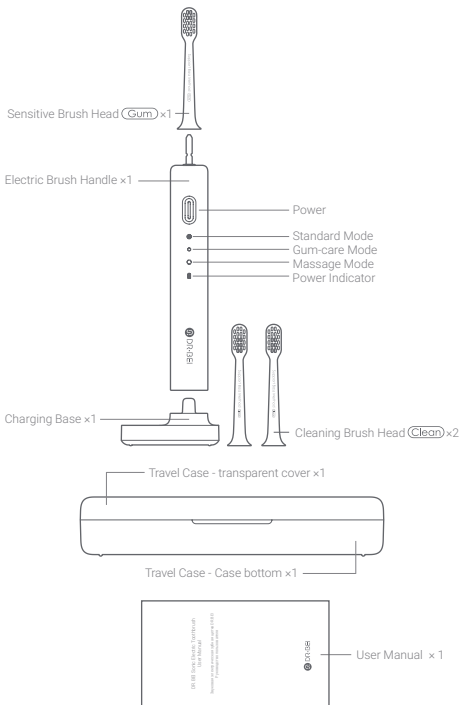
DE.BEI provides you two brush heads for different oral conditions.

Sensitive Brush Head (Gum)

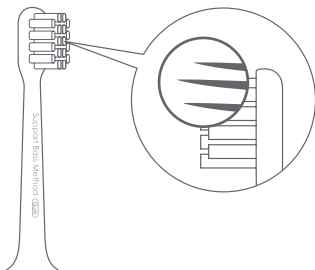Kind reminder *: Sensitive brush head for people with oral sensitivity and gingival bleeding .

Cleaning Brush Head (Clean)

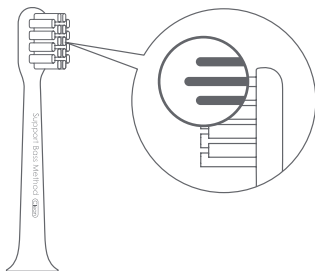Kind reminder *: Cleaning brush head is suitable for people with more teeth loss, such as cigarette stains and tea stains.

*Above instructions are only for your reference. Please choose brush head that brings you the most comfort as oral conditions differs a lot from people to people.

II. How to use

Tip


First time use may accompany itch or slight bleeding due to the fact that your gum just starts getting used to the high frequency cleaning of electric toothbrush. Kindly suggest begin with Gum-care Mode or Sensitive Brush Head and changes later on your needs.

a. Install brush head

Please follow the picture below:

Ensure the bristle point to the heading direction of the button on the handle and then proceed with installation.

Kindly advise that you pick up the handle from the groove beside the power button and remove it from the groove as the bottom of the handle is fixed in travel case.

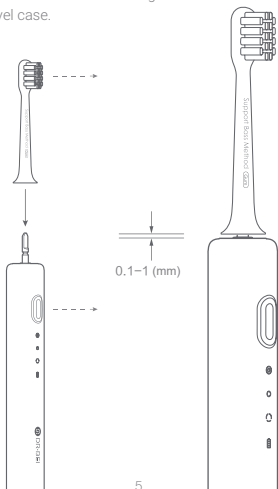
b. Modes Switch

Turn on: Press power button. Default Standard mode for first boot.

Modes Switch: Press power button again to change modes 5 seconds after power on.

Memory mode: The system will record the last mode used for more than 5 seconds as Memory modes and turn to it next time.

Turn off: Press power button to turn off after at least 5 seconds use or automatically turned off after 2 minutes. (3 minutes in massage mode)

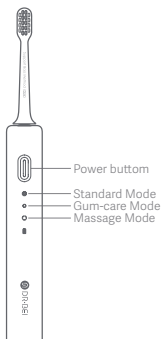c. How to brush Teeth

DR.BEI's Sonic Electric Toothbrush enables a better support for Bass method recommended by American Dental Association (ADA) and effective cleaning of plaque in your mouth.

The following suggestions can help to deliver superior cleaning results :

Bristles and teeth forms a 45° angle can better clean the plaque in gingival sulcus.

Dual frequency and low noise Brushless Magnetic Levitation Sonic Motor


Smart reminder triggered every 30 seconds
Please change to another area after a pause

Gradual increase of vibration
No splashes exerted Improved comfort and cleanness

Kindly advise that you press the brush against your teeth with a force less than 2N and move at every 2 teeth in your first use. Plus with the 30s smart reminder, you will get a thorough cleaning of both in and out side of your teeth. (Note: The massage mode activates the smart reminder every 45 seconds.)

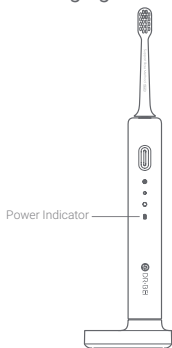d. Charging

 **Power indicator green:** 20%-100% power left, feel free to use.

 **Power indicator orange:** 10%-20% power left, please charge as soon as you can.

 **Power indicator green:** 0%-10% power left, please charge immediately.


While charging, the power indicator is blinking.

Power indicator will stop flash and stay bright when fully charged.

*Cautions:

1. Do not immerse the charging base in water or any liquid.
2. Do not use or place charging base in a high temperature environment.
3. Recommend use of power adapter with CE, UL or local safety standard 5V  1A DC output (need to purchase by the user).
4. Charging protection design adopted. Do not charge on the base for too long in order not to hurt the battery use life .
5. Do not use while charging.


e. Cleaning

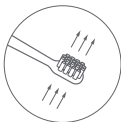Clean the brush head and keep it ventilated

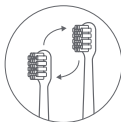Whole body washable
Able to wash under running water

Change brush head every 2-3 months

III. Specifications

| | | |
|--------------------------|------------------|---|
| Electric Brush Handle | Product model | BY-V12 BET-S01 BET-S03 (Violet Gold/Black Gold) |
| | Rated voltage | 3.7V  |
| | Rated power | 2W |
| | Waterproof level | IPX7 |
| | Charging time | At least 4 hours from empty to fully charged |
| Charging Base | Input parameter | 5V  1A |
| Production date | | Printed on the package |
| Production place | | China |

IV. Troubleshooting

| Problems | Reasons | Solutions |
|---|--|--|
| Brush head cannot attach to power shaft correctly | Stuff clog in the brush head socket; Something stick with the power shaft; Brush head is not recommended by manufacturer or not matched; Incorrect orienting. | Remove the stuck matter; Purchase genuine brush head recommended by manufacturer; Align the parts correctly according to the instructions; |
| Cannot turn on | Low power or powered off or in charging. | Charging as soon or remove brush from the charging base. |
| Weak brush head vibration | Incorrect attach or loose connection; Stuff clog between the brush head and power shaft. | Attach correctly; Remove the blocking stuff. |
| Handle cannot connect to charging base | Something stuck at the bottom of handle; Something stuck at the bottom of the socket. | Remove the blocking stuff. |
| Power indicator not showing charge while attach to base | Incorrect connection of plug and outlet. | Check if plug connect outlet correctly. |
| Power indicator cannot show charging state | Low power or used up. | Charging correctly will soon restore it correctly. |

V. After-sales service

If malfunction in DR.BEI Sonic Electronic Toothbrush Performance Failure List occurs, contact local authorized distributor or retailers.

DR.BEI Sonic Electronic Toothbrush Performance Failure List

| Name | Performance Failure |
|--|--|
| DR.BEI Sonic Electronic Toothbrush | Cannot turn on or off |
| | Brush head cannot connect to power shaft when there is no foreign matter |
| | Abnormal vibration of brush head |
| | Power indicator does not flash |
| | Cannot charge the appliance |

WARRANTY EXCLUSIONS

1. Damage caused by unauthorized repair, misuse, collision, neglect, abuse, liquid intake, accidents, alterations or incorrect use of product accessories that is not recommended by the manufacturer, or tear and alter labels or anti-counterfeiting marks.
2. No longer in the validity period of warranty service(repair, refund, replacement).
3. Damage caused by force majeure.
4. Failure to comply with the malfunction listed in DR.BEI Sonic Electronic Toothbrush Performance Failure List.
5. Failure of the product and its accessories in the DR.BEI Sonic Electronic Toothbrush Performance Failure List but caused by users.

VI. Environmental protection

When DR.BEI Sonic Electric Toothbrush reaches its useful life, you may need to remove the internal batteries and recycle them properly according to national or local laws and regulations. If you are unable to remove the battery, we strongly advise you to take your product to an official collection point designated by government.

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Manufacturer: Wuxi Qinghe Xiaobei Technology Co., Ltd.