
DR.BEI Sonic Electric Toothbrush Q3 User Manual

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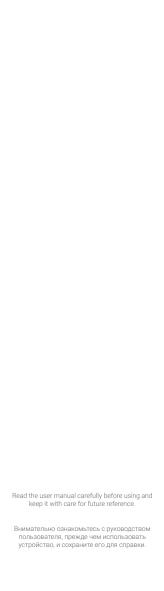Important to Note

△Statement

L'Istatement
First time use may accompany itch or slight bleeding due to the fact that your gum just starts getting used to the high frequency cleaning of electric toothbrush. Kindly suggest begin with Soft Mode. Generally speaking, this phenomenon will disappear after 1–2 weeks. If you do not feel any significant improvement after more than 2 weeks, please stop using this product and promptly consult your dentist. This product is not suitable for users that under the age of 6.

* Tips

- Do not share brush head with others.
- Do not wrap the charging cable the product when it is powered on.

* Notes

- 1. This product is designed for oral cleansing. Do not use it for any other purpose. If you experience any discomfort or pain, please stop using the product and consult your dentist.
- Make sure to check this product for any apparent damage before using it. If damaged or dropped, the product should not be used.
- Before connecting to the power supply, please check whether the local power supply voltage meets the specification of the product.
- Do not use the product if any damage is found on the charging cable or charging port.

* Warnings

- Consult your dental professional before use if you once had orthodontic treatment, dental prosthesis or implantation.
 This product complies with the safety standards for
- electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer before use.

 3. This product contains a built-in rechargeable
- battery. Please keep the product away from flames, and avoid charging, using, or placing it in high temperature.

 4. Make sure the power plug and interface are completely dry before connecting to the power
- supply.

 5. Do not bump or drop the product. If it is damaged due a collision or being dropped, please stop using
- the product.
 6. Do not modify, disassemble or repair the product by
- yourself.
- 7. Do not use the product in shower or bath. 8. Do not use the product if any damage is


. their guardian.

found.(brush head, handle)

9. Anyone without full civil capacity or lack of relevant experience and knowledge (including young children) must not use, clean or maintain the product without the supervision and guidance of

I. Introduction

Thank you for choosing DR.BEI as your oral health care assistant. Please check all the relevant parts are provided as shown in the following list.

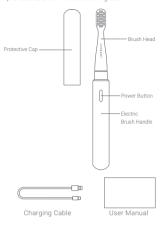
II. How to Use

Please note:

When using this product for the first time, you may experience slight itching or bleeding. This occurs when gums are not accustomed to high-frequency cleaning with an electric toothbrush. It is recommended to use this product in Soft Mode initially and switch modes as required after a period of time.

a. Usage

Remove the protective cap to check that the brush head is well fixed and the gap between the head and the handle is 0.3–0.5 mm.

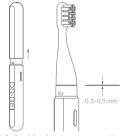Wet the brush head, squeeze an appropriate amount of toothpaste onto the head and begin cleaning.

b. Mode Switch

This product features Cleaning, Soft and Massage modes.

The Cleaning Mode is suitable for daily cleaning, removing plaque and maintaining oral hygiene; the Soft Mode is for those with oral hypersensitivity; and the Massage Mode is for thoroughly removing obstinate tooth stains.

Press the power button to turn th when powering on for the first time. Switch Modes:

Press the power button again within 5 seconds of powering on to switch mo Memory Mode:

When a mode lasts for more than 5 seconds, the system will record this mod as the memory mode. The memory mode will be started by default the next time Power Off-In any situation, press and hold* the

many studation, press and note: the power button to turn the power off. Alternatively, press the power button to shut down 5 seconds after use, or press the power button repeatedly until the power turns off. The brush will also shut down automatically after 2 minutes of use (3 minutes on Massage mode).

ration to press and hold the power button should be more than 1.5 se

c. How to Brush Teeth

Please refer to the following suggestions to improve

better clean the plague in gingival sulcus.

Smart reminder triggered every 30* seconds. hange to another area after a pause.


is recommended to press on teeth with a force below 2N to clean the teeth

When using the product for the first time, it is recommended to press gently on the teeth with a force of less than 2N, cover every area of the mouth, and comprehensively clean both the inside and outside of the teeth. The smart reminder function should be triggered every 30* seconds.

*In Massage Mode, the smart reminder is triggered every 45 seconds.

d. Charging

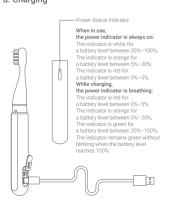*Notes for Charging:

- 1. When charging, please first remove the protective plug of the charging port at the lower part of the product, then insert the small end of the charging cable into the charging port before connecting the other end to a 5V==0.5A DC output power adapter conforming to CE/ETL/KC/PSE, etc. standards (power adapter purchased separately).
- 2. This product is designed with a charging protection function to automatically stop charging when the battery level reaches 100%. Natural consumption of electricity will lead to repeated charging, thus impacting the service life of the battery, so please
- do not leave it charging for a long time. 3. This product cannot be used while charging.

e. Cleaning

Clean the brush head with running water after use.

Change the brush her every 2-3 months.

III. Specifications

Electric Brush Handle	Product Model	DR.BEI Q3
	Rated Voltage	3.7V
	Rated Power	0.6W
	Waterproof Level	IPX7
	Charging Time	At least 4 hours from empty to fully charged
	Rated Input	5V===0.5A
Production Place		China

IV. Troubleshooting

Problems	Reasons	Solutions
Brush head cannot attach to power shaft correctly.	Stuff clog between the brush head and power shaft; Brush head is not genuine or not matched; Incorrect orienting.	Remove the blocking stuff; Purchase genuine brush head recommended by manufacturer; Align the parts correctly according to the instructions.
Cannot power on.	Low power; Product in charging.	Charge the product; Unplug the charger.
Weak brush head vibration.	Incorrect or loose attachment of brush head; Stuff clog between the brush head and power shaft.	Attach correctly, Remove the blocking stuff.
Handle cannot connect to charging cable.	Stuff clog in the charging port.	Remove the blocking stuff.
The brush handle not showing charging status while attaching to the charger.	The charger failed to connect to the power supply.	Check the power supply.
Power indicator does not work.	Low power or used up.	Charging correctly will soon restore it correctly.

V. After-sales Service

If malfunction in DR.BEI Sonic Electronic Toothbrush Performance Failure List occurs, contact local authorized distributors or retailers.

DR.BEI Sonic Electronic Toothbrush Performance

Name	Performance Failure
DR.BEI Sonic Electronic Toothbrush	Cannot turn on or off.
	Brush head cannot connect to power shaft when there is no foreign matter.
	Abnormal vibration of brush head.
	Power indicator does not work.
	Cannot charge the appliance.

WARRANTY EXCLUSIONS


- Damage caused by unauthorized repair, misuse, collision, neglect, abuse, liquid intake, accidents, alterations or incorrect use of product accessories that is not recommended by the manufacturer, or tear and alter labels or anti-counterfeiting marks.
- No longer in the validity period of warranty service (repair, refund, replacement).
- Damage caused by force majeure.
- Fail to comply with the malfunction listed in DR.BEI Sonic Electronic Toothbrush Performance Failure List.
- DR.BEI Sonic Electronic Toothbrush Performance Failure List but caused by users.

VI. Environmental Protection

When the product reaches its usage expectancy, according to the laws and regulations of the country or region, it may be necessary to remove the internal battery in order for it to be recycled and properly handled. If you are unable to remove the battery, please send the product to a government-designated recycling center.

Instructions and suggestions regarding battery waste disposal:

It is not recommended that the battery is replaced by any user. Forcing the battery out will damage the product's functionality. Please exercise with caution.

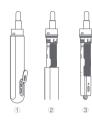- ① Open the protective cap of the USB port and remove the screw.
- Remove the shell from the handle.
 Cut the nickel tab behind the battery and remove the battery.

The illustrations of the product, accessories, user interface, etc., in the manual are schematic diagrams for reference only. Due to the updating and upgrading of products, there may be slight differences between the actual product and the schematic diagram. The actual product may be a more up-to-date version.

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