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IMPORTANT SAFEGUARDS

⚠ Statement

While using the product for the first time, you might experience a tickling, tingling, or slight bleeding as your gums are not used to such efficient cleaning. You are advised to start from the soft mode. Normally, the sensation subsides after 1 or 2 weeks of the use; if not, you should immediately stop using the product and consult your dentist immediately.

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING

* **DANGER**– To reduce the risk of electrocution:

1. Always unplug this product immediately after using;
2. Do not use while bathing;
3. Do not place or store product where it can fall or be pulled into a tub or sink;
4. Do not place in or drop into water or other liquid;
5. Do not reach for a product that has fallen into water. Unplug immediately.

* **WARNING**– To reduce the risk of burns, electrocution, fire, or injury to persons:

1. A product should never be left unattended when plugged in;
2. Close supervision is necessary when this product is used by, on, or near children or invalids;
3. The product is designed only to clean the mouth. Do not use it for any other purpose. In case of pain or any other complaint, please stop using the product and consult your dentist;
4. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the product to a service center for examination and repair;
5. Keep the cord away from heated surfaces;
6. Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and the like;
7. Never use while sleeping or drowsy;
8. Never drop or insert any object into any opening or hose;
9. Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered;
10. Unplug this product before filling. Fill (reservoir) with water only unless otherwise specified by manufacturer. Do not overfill;
11. Before connecting to power, make sure that supply voltage is matched with the rated voltage of the product;
12. Before maintenance, make sure that the product is not being charged and is powered off;
13. Do not fill undiluted mouthwash into the product; mouthwash should be diluted according to the relevant specifications before filling it into the product;
14. People undergoing orthodontic treatment or dental filling, or have dental implants should use the product only after consulting a dentist;
15. Do not keep the product in a hot or damp place; keep it away from direct sunlight;
16. The product has an in-built rechargeable battery. Do not expose it to flames or heat it, or charge, use, or place it at high temperatures;
17. Make sure that the power plug and the port are completely dry before putting the product on charge;
18. Do not try to modify, disassemble or repair the product on your own;
19. Do not fill the water tank with water above 40°C (104°F) temperature for cleaning your mouth as it may cause injuries;
20. Do not press the nozzle hard against the teeth or gums as it may cause injury;


* **Tips**

1. Do not share nozzle with others;
2. Do not wrap the cord around the product when it is switched on;
3. Do not use the product without nozzle installed;
4. Household use only

Instructions and recommendations for the handling of a used battery:

When the product has reached the end of its service life, it may be necessary to remove the internal battery for recycling according to the prevailing national or local laws and regulations. If you cannot take out the battery, please hand over the product to the recycling center as designated by the government.

The battery of the product should not be replaced by users on their own; forcible removal may result in damage to the product, please proceed with caution.

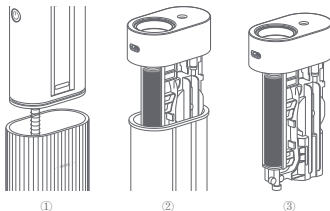
- ① Remove the water tank;
- ② Remove the lower shell of the product;
- ③ Cut off the battery nickel bar behind the support and remove the battery.


The illustrations of the product, parts and user interfaces in the instruction manual are for information only. The products are subject to change and upgrade, and thus may differ slightly from the illustrations.

I. Product Description

Thank you for choosing DR.BEI as your oral health assistant. Please check if the product is complete with relevant parts mentioned in the list below.

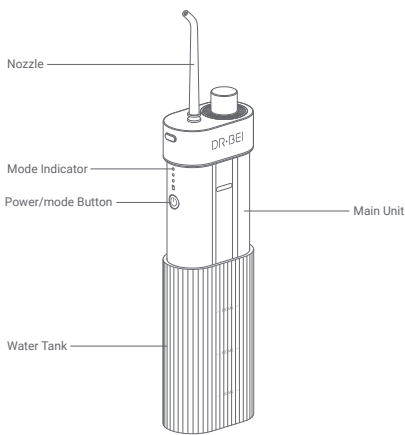User Manual

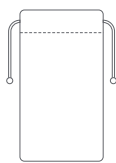Travel Bag


USB Cable

II. Usage


Tips
Before use, check if the product has sufficient power. Before first use, charge the product. A red light indicates that the battery power is low and the product should be charged.

a. Install or remove the nozzle

1. Open the flip top and take the nozzle out of the storage compartment.

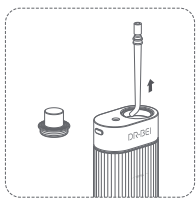
2. Insert the nozzle firmly into the nozzle socket until it clicks into place.

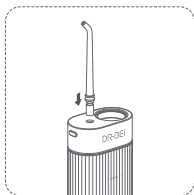
3. Press and hold the nozzle release button, pull out the nozzle gently, and restore it into the storage compartment.

Warning: Do not press the nozzle release button when the product is running.

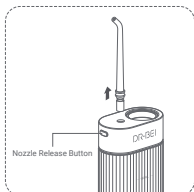
b. Fill with Water


1. Open the flip top and pull down the water tank part to the usage line.

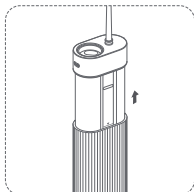2. Keep the product upright and pour water in it.

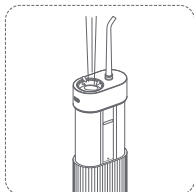
3. Close the flip top once adequate quantity of water is filled.

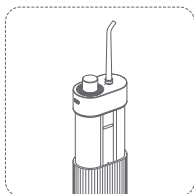Clean tap water or purified water is preferred; lukewarm water not over 40 °C (104 °F), saline solution with a salt concentration of less than 5% or diluted mouthwash are also acceptable.

c. Switch Between Modes

The product has a Cleaning Mode, a Soft Mode, and a Massage Mode.
The Cleaning Mode is for normal cleaning and applicable to most people; the Soft Mode is for new users of water flossers or those with sensitive teeth; with the Massage Mode, teeth are cleaned and gums are massaged.

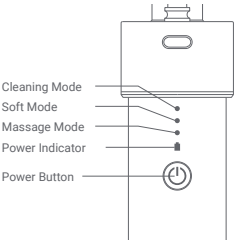Power On:
Press down the power button to switch on the product. For first use, the product will work in the cleaning mode by default.

Switch between modes:
To switch between modes, press the power button again in 10 seconds after the product is switched on.


Memory Mode:
When certain mode is used for over 10 seconds, the product will store the mode as the memory mode. It will start in the same mode at the next boot.

Shut Down:
To switch off the product, press and hold the power button.
Alternatively, press the power button to switch off the product after it has run for 10 seconds;
or long press the power button until the product shuts off;
or the product shuts off automatically after 2 minutes.

* Press the power button and hold for at least 1.5 seconds.

d. Clean Teeth

1. Place the nozzle in your mouth and aim at your teeth, place it between teeth or gums; make sure that the water flow is perpendicular to the part for flushing and your mouth is slightly open.

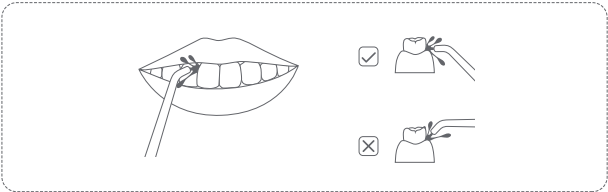
2. Press down the power button for cleaning the teeth. During the use, keep the product upright to prevent overflow of water as the water suction tube inside the water tank may come out from the water.
3. The product will shut off automatically after it has run for 2 minutes; if you wish to use it for some more time, press the power button.


e. Clean and maintain

Clean After Using the Product
There is a risk of bacteria propagating, so you are advised to wash and dry the product thoroughly before storing it.


1. Switch off the product. Open the flip top and dispose of the water in the water tank.

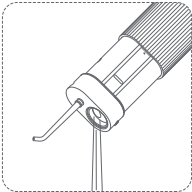2. Press the power button again to switch on the product and drain the water inside the water tank.

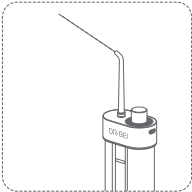3. Push up the water tank to the top to restore the main unit; place the nozzle back into the storage compartment; close the flip top.

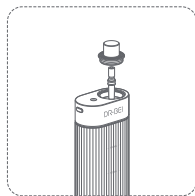Maintain after using the product

Nozzle

Wash with water and wipe with a clean and dry cloth. It is recommended to replace the nozzle once in 6 months; if the nozzle gets deformed, please replace it with a new one; do not bend, stretch or twist the nozzle in any case.

Main unit

The main unit can be wiped with a soft dry cloth or a damp cloth. Do not immerse the product in water for washing.

Water suction tube


Wash with water and wipe with a cloth. Do not bend the water suction tube.

Filter


Wash with water or soft brush and dry with a cloth. Do not remove it while washing.

Water tank

The water tank is removable for easy cleaning. It is recommended to clean the water tank weekly.

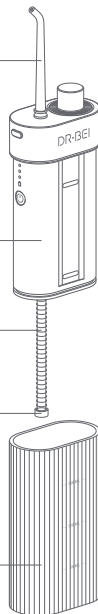f. Charge

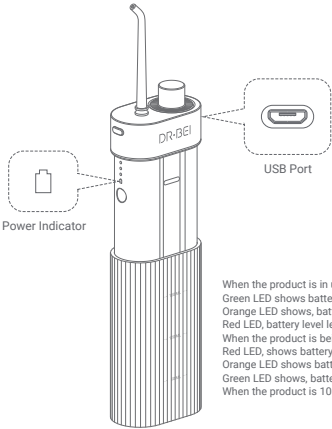When the product is in use, the power indicator is on:
Green LED shows battery level between 20%-100%;
Orange LED shows, battery level between 10%-20%;
Red LED, battery level less than 10%;
When the product is being charged, the power indicator is winking:
Red LED, shows battery level less than 10%;
Orange LED shows battery level between 10%-20%;
Green LED shows, battery level between 20%-100%;
When the product is 100% charged, the indicator is in solid green.

* CAUTIONS FOR CHARGING:

1. Before charging, plug the smaller end of the cord into the USB port, and plug the other end into an CE/ETL/KC/PSE, etc. approved 5V DC 1A adapter (to be purchased separately);
2. The product is designed with charging protection where charging gets stopped when the battery is 100% charged. Please do not keep charging the battery for longer time as it will affect the lifespan of the battery adversely due to natural drain and repeated charging;
3. The product cannot be used when it is charging.

III. Specifications

Product Name		DR.BEI Portable Water Flosser
Water flosser	Model	DR.BEI GF3
	Rated Voltage	3.7V DC
	Rated Power	5W
	IP Rating	IPX7
	Charging Time	Full recharge takes at least 4 hours
	Rated Input	5V DC 1A
Place of Origin		China

IV. Troubleshoot

Malfunctions	Causes	Solutions
The product does not operate even when the switch is turned on.	The battery power is low, or the battery is empty or being charged.	Charge the battery or remove the product from the charger.
The water pressure is low.	The nozzle is clogged or damaged.	Replace the nozzle.
The indicator is abnormal or does not light up.	There is a circuit fault.	Return the product to the after-sales service center.
The product does not indicate charging even when the cord is connected.	The plug is not properly inserted into the electric outlet.	Check if the electric outlet is connected properly.

V. After-sales Service

If the product has any fault mentioned in the List of Product Faults, please contact a local dealer or an after-sales service provider.

Warranty restrictions (the warranty does not cover the following)

- 1. Unauthorized repair, misuse, damage by shock, neglect, abuse, permeation with liquid, accident, or alteration, or incorrect use of any part not for the product, or torn or altered labels or anti-counterfeit marks;
- 2. Expiry of warranty;
- 3. Damage due to force majeure;
- 4. Inconsistency with the faults mentioned in the List of Product Faults;
- 5. Faults of the product or its part which is resulted from deliberate damage and not mentioned in the List of Product Faults.

List of Product Faults

Name	Faults
DR.BEI Portable Water Flosser	The product fails to switch on or off.
	No water comes out when the water tank has water and there is no dirt stuck in the nozzle.
	The power indicator is not functioning when the battery is being charged.
	The battery does not charge properly.

After-sales Service: service@doctorbei.com
Manufacturer: Wuxi Qinghe Xiaobei Technology Co., Ltd.
Address: Room 1106, 11th Floor, NO.2 of Building 530, Binhu District, Wuxi City, Jiangsu, P.R. China

VI. Environmental Protection

When the product has reached the end of its service life, it may be necessary to remove the internal battery for recycling according to relevant national or local laws and regulations. If you cannot take out the battery, please hand over the product to the recycling center designated by the government.

EU Declaration of Conformity

We Wuxi Qinghe Xiaobei Technology Co., Ltd., hereby, declares that this equipment is in compliance with the applicable Directives and European Norms, and amendments.

WEEE Disposal and Recycling

Correct Disposal of this product. This marking indicates that this product should not be disposal with other household wastes throughout the EU.

To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources.

To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

Warning: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/ TV technician for help.

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For more information about Dr. BEI's products,
please visit DR.BEI's official website:

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Manufacturer: Wuxi Qinghe Xiaobei Technology Co., Ltd.
Address: Room 1106, 11th Floor, NO.2 of Building 530, Binhu District, Wuxi City, Jiangsu, P.R. China
After-sales service: service@doctorbei.com

Fabricante: Wuxi Qinghe Xiaobei Technology Co., Ltd.
Dirección: Room 1106, 11th Floor, NO.2 of Building 530, Binhu District, Wuxi City, Jiangsu, P.R. China
Servicio postventa: service@doctorbei.com

Производители: Wuxi Qinghe Xiaobei Technology Co., Ltd.
Адрес: Room 1106, 11th Floor, NO.2 of Building 530, Binhu District, Wuxi City, Jiangsu, P.R. China
почтовый ящик: service@doctorbei.com

Hersteller: Wuxi Qinghe Xiaobei Technology Co., Ltd.
Adresse: Room 1106, 11th Floor, NO.2 of Building 530, Binhu District, Wuxi City, Jiangsu, P.R. China
Kundendienst: service@doctorbei.com

Produttore: Wuxi Qinghe Xiaobei Technology Co., Ltd.
Indirizzo: Room 1106, 11th Floor, NO.2 of Building 530, Binhu District, Wuxi City, Jiangsu, P.R. China
Assistenza post-vendita: service@doctorbei.com

Fabricant: Wuxi Qinghe Xiaobei Technology Co., Ltd.
Adresse: Room 1106, 11th Floor, NO.2 of Building 530, Binhu District, Wuxi City, Jiangsu, P.R. China
Assistenza post-vendita: service@doctorbei.com

メーカー：Wuxi Qinghe Xiaobei Technology Co., Ltd.
アドレス：Room 1106, 11th Floor, NO.2 of Building 530, Binhu District, Wuxi City, Jiangsu, P.R. China
カスタマーサービス：service@doctorbei.com

Scan to connect and
discover more on Facebook.

 DR.BEI